


**Variable Speed Vortex Mixer  
MS-VM Series**

↓ ↓ ↓  
**Auto Off Touch**



MS-VM-001

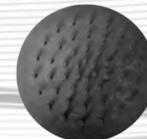

MS-VM-002



MS-VM-003



MS-VM-004



MS-VM-005




MS-VM-006



MS-VM-007



MS-VM-008

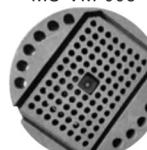

MS-VM-009



MS-VM-010



MS-VM-011



Model	Variable Speed Vortex Mixer
Cat. No.	MS-VM Series
Voltage	AC 110V or 220V, 50/60 Hz
Power	60W
Movement	Orbital
Orbital Diameter	4mm
Speed Range	0-2800rpm, manual adjustment
Mode	Continuous(ON), Depressing(Touch), Stop(OFF)
Dimensions	W170×D120×H170mm
Weight	4.8kg
feet	4 Rubber Feet Pad
Standard Accessory	3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head
Optional Accessory	15ml/50ml Tube Holder
	1.5ml/15ml/50ml Tube Foam Rack
	96-well PCR Plate Foam Platform
	96-well Microplate Foam Platform

The major science Vortex Mixer series is an efficient and compact instrument for any routine mixing applications needed in the field of life sciences, chemistry, medicine, bioengineering, analytical and research laboratories.

This product includes the eccentric rotating structure to generate vortex flow in the liquid in test tubes centrifuge tubes, microplate and other containers and offers fully mixing the solution.

Choose from two modes of operation : continuous mode provide variable speed can allows low rpm start-up for gentle mixing or high speed mixing for vigorous vortexing of samples, touch mode which activates mixing when depressing the platform.

The major science Vortex Mixer series is stable and reliable with reliable durability and also has a variety of accessories to choose from to adapt to different experimental needs.

Cat. No	Product Description
MS-VM-110-US	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 110V, US Power Cord
MS-VM-110-JP	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 110V, JP Power Cord
MS-VM-220-EU	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 220V, EU Power Cord
MS-VM-220-UK	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 220V, UK Power Cord
MS-VM-220-IN	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 220V, IN Power Cord
MS-VM-220-CN	Variable Speed Vortex Mixer with 3" Round Rubber Flat Cover & Pop-Off Rubber Cup Head, 220V, CN Power Cord

Cat. No	Product Description
MS-VM-001	3" Round Rubber Flat Cover
MS-VM-002	Pop-Off Rubber Cup Head
MS-VM-003	15ml Tube Horizontal Holder, 4 clamps
MS-VM-004	50ml Tube Horizontal Holder, 2 clamps
MS-VM-005	1.5ml Eppendorf / Tube Foam Rack
MS-VM-006	15ml Tube Foam Rack
MS-VM-007	50ml Tube Foam Rack
MS-VM-008	96-well PCR Plate Foam Platform
MS-VM-009	96-well Microplate Foam Platform
MS-VM-010	15ml Tube Horizontal Holder, 2 clamps
MS-VM-011	50ml Tube Horizontal Holder, 4 clamps



\* All images are for reference only, actual products might differ from the pictures above.  
\* Technical specifications subject to change without notice.